


■ Made
■ in
■ Germany

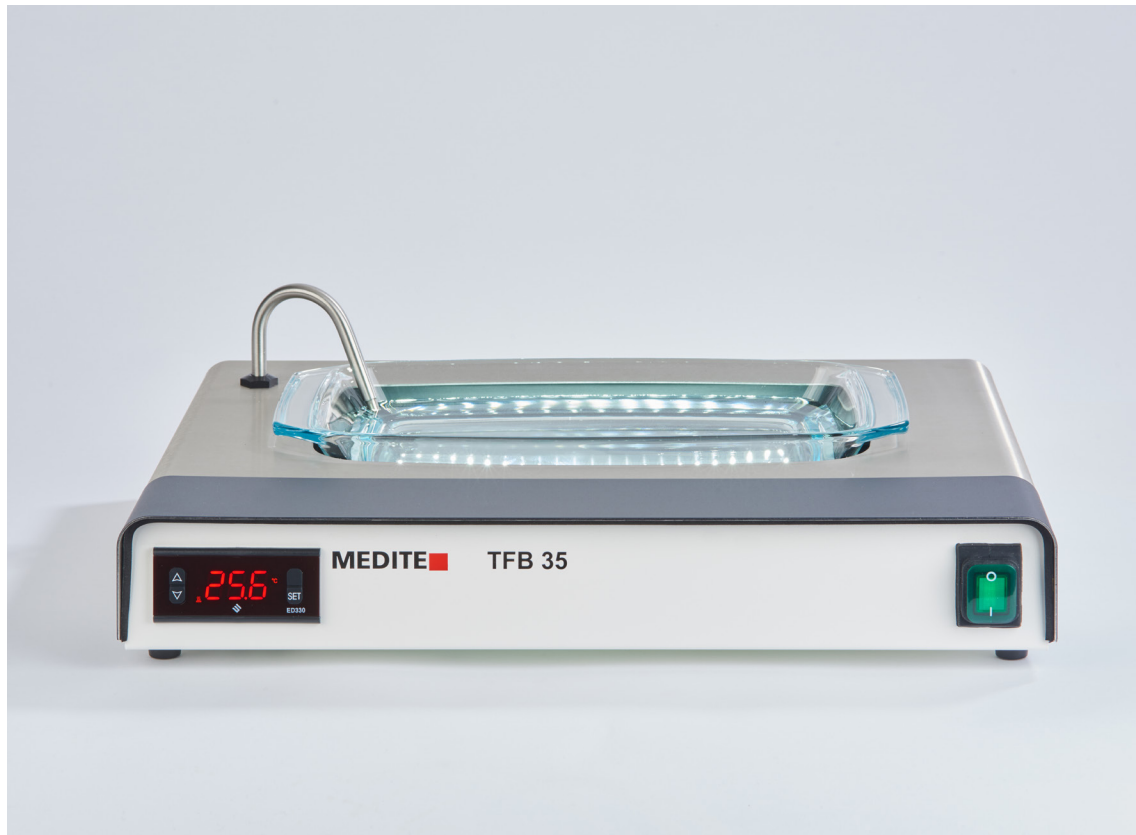

The Tissue Floatation Bath **TFB 35** is an LED-illuminated, rectangular stretching bath for paraffin sections that has been tried and tested in the routine of everyday laboratory work. The central component is the removable glass tray, under which a LED light strip is located. The light comes from the rear edge under the glass tray, ensuring excellent visibility of the paraffin sections on the water surface.

The glass tray of the **TFB 35.000** has a volume of 1.5 liters that can be easily removed and also cleaned in the dishwasher. A second model, the **TFB 35.500**, has a higher glass tray with a capacity of 2.8 liters.

In addition to a removable glass tray, the **TFB 35** has a heated work surface on the front. It is coated with a resistant black plastic and offers enough space to store slides for drying.

The temperature is set and controlled by a digital controller. To ensure optimal temperature control, the **TFB 35** has a temperature sensor that is immersed in the glass bowl and constantly measures the water temperature.

A matching lid, which ensures temperature stability and protection against dirt when not in use, is available separately.



Technical Data TFB 35.000		Cat. No.
Bath Capacity:	max. 1.5 liters	
Temperature:	+30 °C bis +55 °C	
Dimensions Heated Front Area (W/D):	370 x 75 mm	
Dimensions Glass Tray (W/D/H):	270 x 165 x 50 mm	
Dimensions Overall (W/D/H):	430 x 365 x 160 mm	
Weight:	ca. 7 kg	
Power Supply:	230 V / 50 – 60 Hz / 250 VA	01-3500-02
	115 V / 50 – 60 Hz / 250 VA	01-3510-02
Technical Data TFB 35.500		
Bath Capacity:	max. 2.8 liters	
Temperature:	+30 °C to +55 °C	
Dimensions Heated Front Area (W/D):	370 x 75 mm	
Dimensions Glass Tray (W/D/H):	270 x 165 x 80 mm	
Dimensions Overall (W/D/H):	430 x 365 x 190 mm	
Weight:	ca. 7 kg	
Power Supply:	230 V / 50 – 60 Hz / 250 VA	01-3550-02
	115 V / 50 – 60 Hz / 250 VA	01-3560-02

Accessories	Purchasing Unit	Cat. No.
Lid made of Aluminum for TFB 35.000	1 pc.	11-4550-01
Lid made of Aluminum for TFB 35.500	1 pc.	11-4550-10
Glass Tray, Capacity 1.5 liters, for TFB 35.000	1 pc.	51-4552-00
Glass Tray, Capacity 2.8 liters, for TFB 35.500	1 pc.	51-4553-00

MEDITE

MEDITE Medical GmbH

Wollenweberstr. 12 • 31303 Burgdorf • Germany
 Phone +49 5136 88 84-0 • Fax +49 5136 88 84-55
 info@medite.de • www.medite.de

Subject to change without prior notice. Pictures and colours may vary.



Distributed by:

Axlab A/S

+45 35 43 18 81
 axlab@axlab.dk
 www.axlab.dk

Axlab

YOUR VISION.
 OUR INSPIRATION